

Fitting your hearing aid – Slim Tube

- 1 Place the hearing aid behind your ear with the soft tip facing towards your head.



- 2 Gently push the soft tip into the ear canal with the tail facing backwards.



- 3 Fold the tail backwards into the bowl of the ear.




- 4 Correctly fitted aid.





Please note if you are fitting a right hearing aid use your right hand, and if fitting a left hearing aid use your left hand.

Cleaning Your Hearing Aid

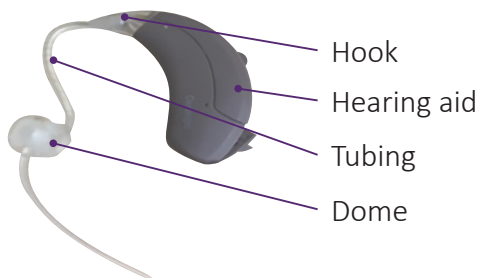

- 1 Unscrew the tubing clockwise from IN FRONT of the hook on the hearing aid. The hook contains a filter which you cannot push the cleaning wire through.
- 2 Gently push the cleaning wire through the tubing from top to bottom. Pull the cleaning wire all the way through.
- 3 Remove any debris from the end of the tubing and wipe the dome to remove any excess wax.
- 4 Gently push the tubing back on to the hook then turn it anti-clockwise until the tubing is back to its original position. Be careful not to kink or over tighten the tube.

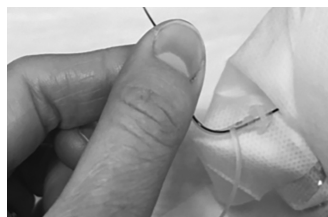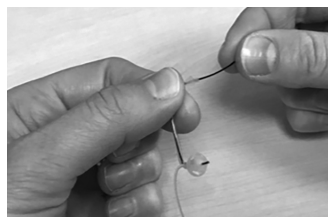

Please note: Tubing should be replaced every 4-6 months